- 2 SEECODE DESIGN-SMART-KÖRPERWAAGE
- 12 SEECODE DESIGN-SMART-BODY SCALE
- 22 SEECODE DESIGN-SMART-PÈSE-CORPS
- 32 SEECODE DESIGN-SMART-LICHAAMSWEEGSCHAAL

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[www.mobiset.de](http://www.mobiset.de)

## Instruction manual

Thank you for purchasing our SEECODE Design Smart Body Scale, which can provide you with a lot of important data about your body. With this scale, you can conveniently measure such things as body weight, heart rate, body balance, BMI, body fat percentage, body weight without fat, subcutaneous fat, visceral fat, body water, bone mass, muscle mass, protein balance, basal metabolic rate (BMR) and biological age.

In order to enjoy your product for a long time and to use it to its full potential, please read these operating instructions carefully before using the scale to avoid damage due to incorrect operation. Please pay particular attention to the SAFETY INSTRUCTIONS. If you give the machine to a third party, this instruction manual must also be handed over.

Failure to follow the instructions in this manual or use of the scale in a manner other than that specified in this manual may result in personal injury or damage to the scale.

## User manual

Please keep this manual for reference. The Design Smart Body Scale is intended for use by adults. It is not a toy and should not be operated by children. Please keep this product out of the reach of children for their safety.

## Package contents

SEECODE Design Smart Body Scale  
Instruction manual

## QR code for Fitdays APP

Scan the corresponding QR code to download the APP. For detailed information, please refer to the rest of the instruction manual.


## Fitdays APP functions

- Weight
- Heart rate
- Body Balance
- BMI
- Body fat percentage
- Fat-free body weight
- Subcutaneous fat
- Visceral fat
- Body water - Skeletal muscle percentage
- Bone mass
- Muscle mass
- Protein balance
- Basal Metabolic Rate (BMR)
- Biological age
- Body growth of adolescents (age < 16 years)

## Display and functions

1. Battery charge level
2. Bluetooth connection
3. Weight unit: kilogram / pound
4. Body weight
5. Body balance
6. Body fat percentage
7. BMI / Body Measure Index
8. Body water percentage
9. Muscle percentage
10. Heart rate
11. Bone mass
12. Compartment for 4 AAA alkaline batteries, LR03
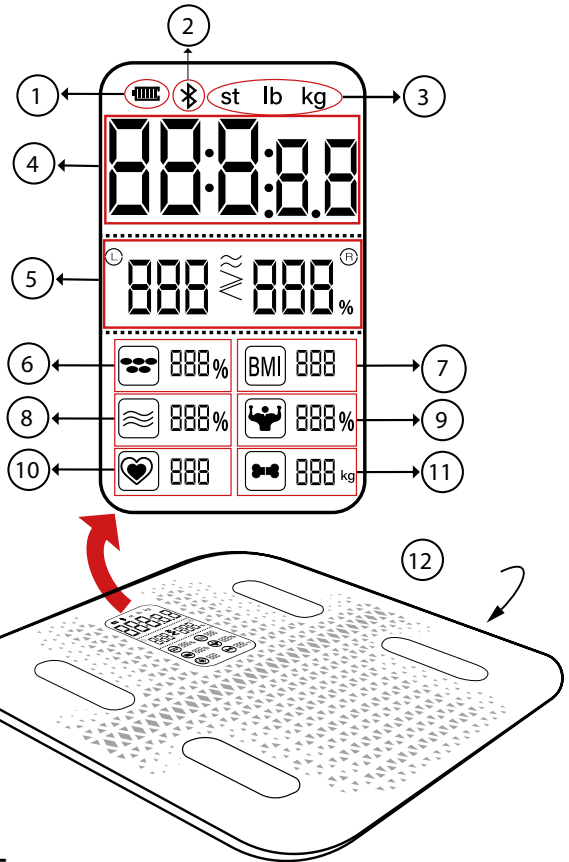
## Notes on the display

-  Battery status indicator
- LO Low battery power
- Err Overload
- C During measurement
- An error occurred
- Step off the scale and back on again to repeat the measurement process.

## SAFETY INSTRUCTIONS, PLEASE NOTE

- Do not use with medical implants, e.g. pacemakers or other electronic medical devices.
- Not recommended for pregnant women. Any measurement results obtained with this device are for reference only and should not be considered as a medical opinion.
- Consult your GP or doctor before making any changes to your diet, exercise plan or other physical activity. Use not recommended for infants, toddlers and children under 10 years of age without adult supervision. Do not stand on the edge of the scale (risk of tipping over) and do not jump on the scale. Do not overload the scale (max. 180kg).
- Do not drop the scale or let any objects fall on it as this may damage the sensors.
- Keep the Design Smart Body Scale away from heat sources such as radiators, etc.
- Do not immerse the scale in water or use chemical cleaning agents.

cleaning agents. Only clean the scale with a slightly damp cloth.  
Make sure that your feet and the scale are dry before stepping on the scale.  
Always place the scale on a firm, dry and flat surface before measuring.  
Do not disassemble, rebuild or repair the individual components of this scale.

Please check for damage to the scale at regular intervals.

## **SPECIFICATIONS / TECHNICAL DATA**

Dimension:	300 x 300 x 25 mm
Display:	VA Screen Display 62 x 109mm
Weight:	1,8 kg
Weight unit:	Kg/Lb
Weight limit:	6 kg - 180 kg
Power supply:	4 AAA alkaline batteries, LR03 (not included)
Number of users:	No limit (each account set up can add 24 people).
Baby mode:	Allows you to weigh the baby's weight
Body growth:	Allows you to record the height of adolescents up to 16 years
Operating temperature:	10 - 40 degrees Celsius
Ambient humidity:	20% - 90%

### **Inserting the batteries**

The 4 AAA alkaline batteries are inserted at the bottom of the scale.

Open the compartment by sliding the cover and insert the batteries. Please make sure that the polarity of the plus and minus poles is correct.

Do not touch any live parts with wet cloths or wet hands.

### **Supported devices for the app:**

iOS 8.0 or higher

Android 6.0 or higher

### **Install the Fitdays app**

1. search for „Fitdays“ in the Apple Store or on Google Play or scan the QR code.
2. download and install the app on your device.
3. the app icon appears on your smartphone or tablet, after the installation is complete.

## Using the Fitdays app


1. place the scale on a firm, level surface.


2. turn on Bluetooth on your smartphone or tablet.
3. register your own ID in the app with your email and a unique password.

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4. add your personal data and then click on ‚Confirm‘. Note: Entering incorrect data such as gender, age, height etc. may result in the measurement being may result in an inaccurate measurement for your person.

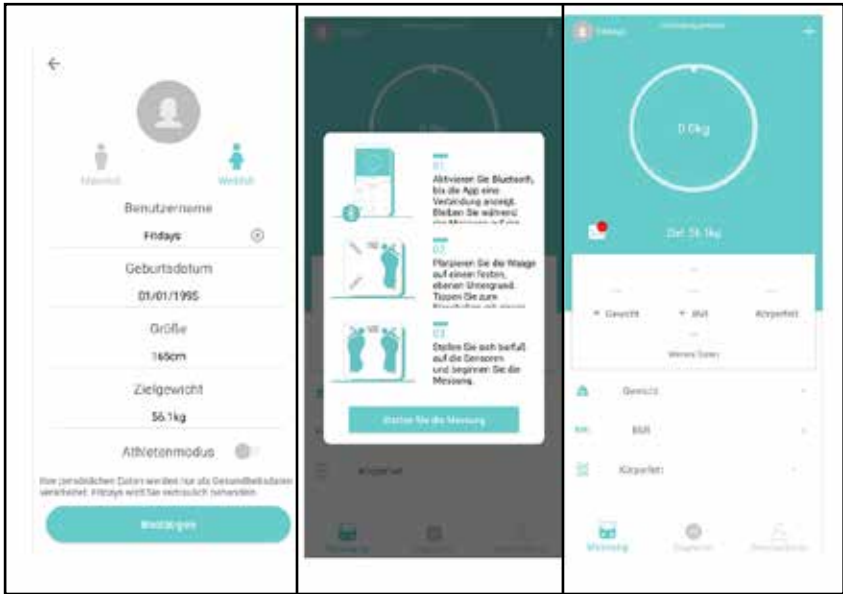- pair the scale via Bluetooth. Click on „User account“ then „Device“, then „+“ to find and pair the scale ID. Note: If the connection is successful, a Bluetooth icon is displayed on the scale and the word will be displayed on the scale and the word ‚Connected‘ will be displayed at the top of the Home screen.

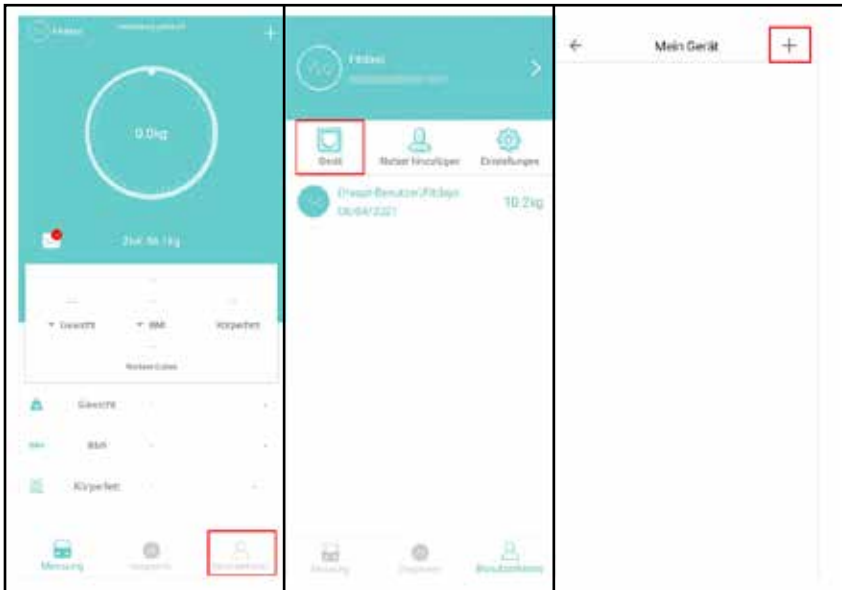6. tap the scale with your foot to activate the display. Wait until the display shows „0,0 „kg.
7. stay at ‚Fitdays‘ on the home screen and step on the scale with bare feet to start the measurement.


Tip:


Stay on the app page, only then can the data be synchronised with the phone.

If you step on the scale with socks, only the weight and body balance will be measured.


8. the number on the display flashes when the scale weighs and analyses you. When the scale shows your left and right weights on the display, you have 5 - 10 seconds to balance your body. Please try to balance your left and right weight equally on the scale.
9. when heart rate is measured, a heart symbol will light up on the scale and flash in the app. Please remain calm and composed to get accurate data.
10. remain on the scale until „Done“ is displayed, then all data on the scale and smartphone will be synchronised. on the scale and the smartphone will be synchronised.

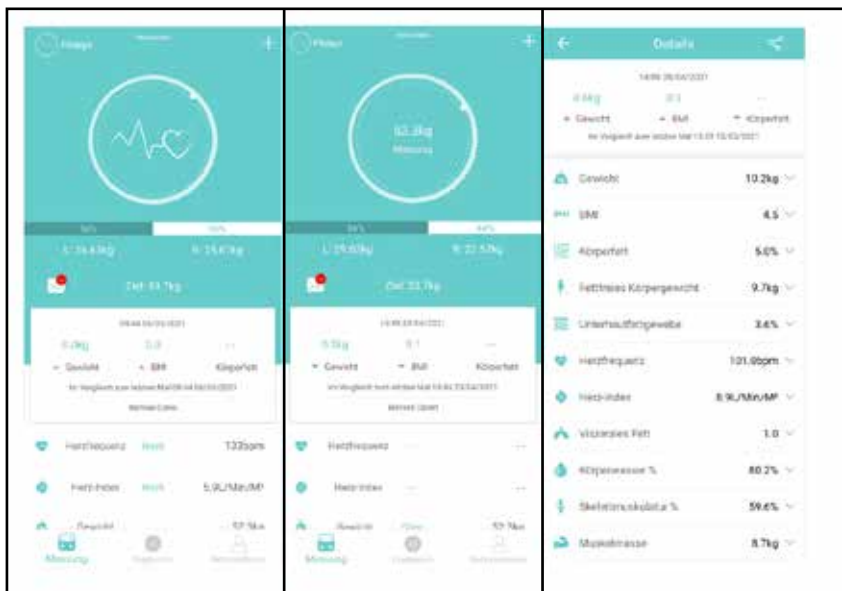
11. click „More data“ to see details of body weight, BMI, BMR, muscle mass, bone mass, etc. check the data difference in the chart by week, month and year in the vertical or horizontal version, control with the icon on the top right.

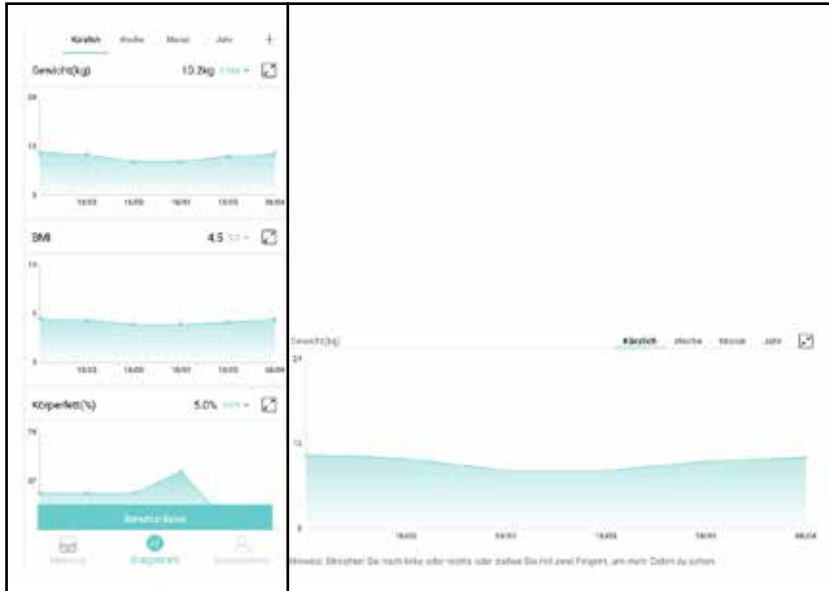12. add or change a user: click on „your profile“ on the homepage or on „User account“, then „Add user“. The smart scale can recognise the user by the weight difference and remind you to select the correct user. Each account can support 24 users

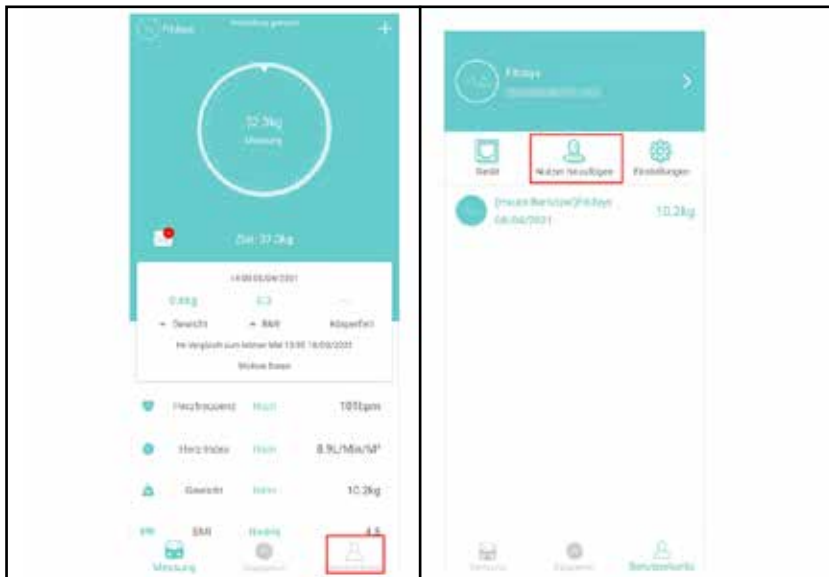13. baby weight check: Click „+ --> baby mode“. An adult steps on the scale first to weigh themselves. After the scale has determined the weight, please step off the scale, hold the baby in your arms and weigh yourself. The Fitdays app will then automatically calculate the baby's weight.
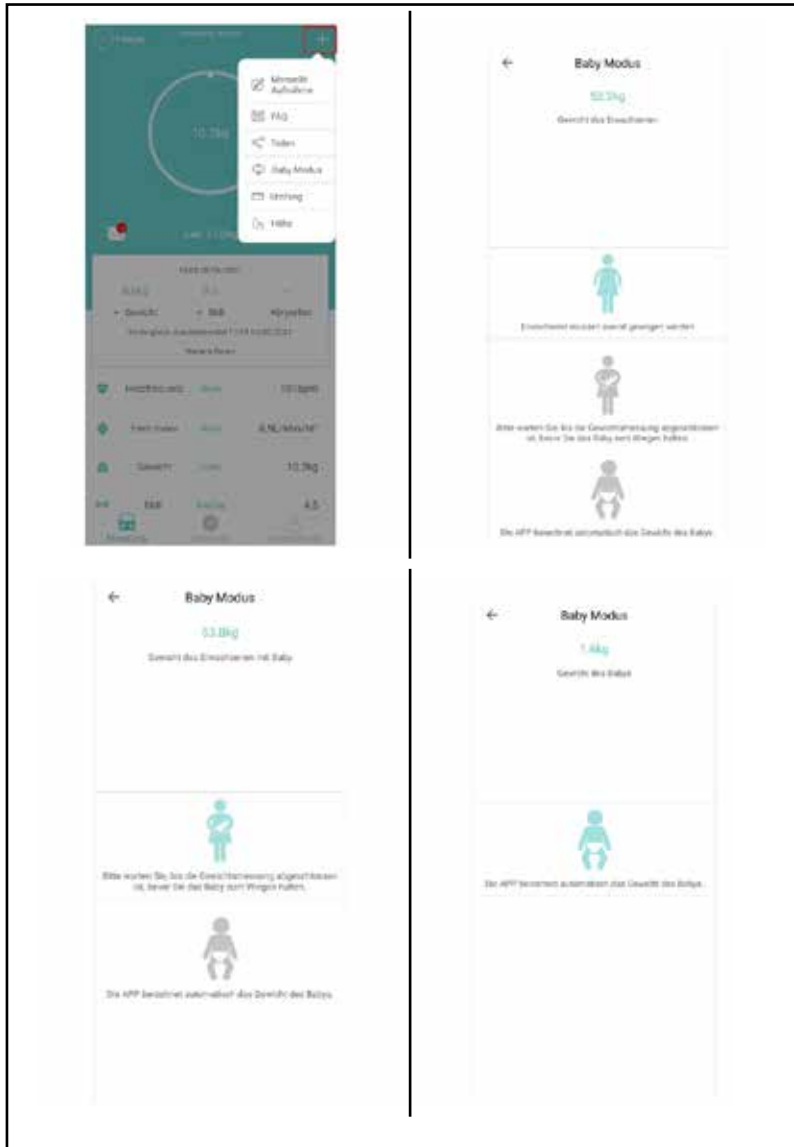
## Troubleshooting guide

### 1. Bluetooth connection failed

Make sure the scale and Bluetooth are turned on on your smartphone.

Check the OS version of your smartphone, it must be iOS 8.0 or higher/Android 6.0 or higher.

If the connection is successful, a Bluetooth symbol will appear on the display and the word.

„Connected“ is displayed in the circle of the smartphone's start page.

If the connection via Bluetooth fails, please click on „My--->Device“ to find the to find the scale ID. 2.

### 2. no body fat data is displayed after weighing.

The user must step on the scale with bare feet.

Make sure that your feet and the scale are dry.

Do not step off the scale if the measurement is not completely finished.

Stand still on the scale for approx. 10 seconds until the number shown on the display stops flashing.

Bluetooth is not turned on or paired with another scale. 3.

### 3. incorrect weight data appears

Check that the scale is on a hard, level surface.

Check each sensor foot on the bottom of the scale so that nothing is stuck to a foot.

### 4. incorrect fat to muscle ratio is displayed in the app.

The body parameter settings are incorrect, check that you have entered the correct gender, height and age.

Check that you have entered the correct gender, height and age.

### 5. recalibration of the scale after transport or placement in another location.

Step on the scale to switch it on. Wait a while and let the scale automatically reset to 0.0 kg. 6.

### 6. The scale will not turn on.

Check if the batteries are empty and insert new ones if necessary.

### 7. exchange another scale on the mobile device.

Click on „User account“, then „Settings“ and „Logout“ to then connect to a new scale.

## Important recycling information

Please send the device to certified drop-off points for recycling.

Before transporting the unit, make sure that it is switched off and use the original packaging if possible.

Used batteries must not be disposed of in household waste because they may contain toxic

components and heavy metals that can be harmful to the environment and health. Give empty batteries to a suitable recycling facility.

### **Guarantee/warranty**

The warranty period is 24 months from the date of purchase. This includes 6 months warranty from the date of purchase. Proof of purchase is the sales receipt. The product has been thoroughly checked for quality and functionality before dispatch.

The warranty is excluded in the following cases:

in the event of faulty installation

in the event of modifications to the system without our approval

damage caused by moisture on the live parts, damage caused by dropping or damaged cables if there is no proof of purchase.

cables if no proof of purchase is available

### **CE Declaration of Conformity**

Mobiset GmbH hereby declares that this product is in conformity with the essential requirements of the applicable European directives.

### **General information**

#### **Copyright**

This document is protected by copyright. Any duplication or reprinting, including excerpts, as well as the reproduction of illustrations, even in modified form, is only permitted with the written consent of the manufacturer.

#### **Notes on environmental protection**

The packaging materials used can be recycled. Dispose of packaging materials that are no longer needed in accordance with the locally applicable regulations.

This product must not be disposed of in normal household waste within the European Union. Dispose of the appliance at municipal collection points.

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Mobiset GmbH – Rösrather Straße 333 – 51107 Köln