
Seecode Thermo Insektenstich-Gerät

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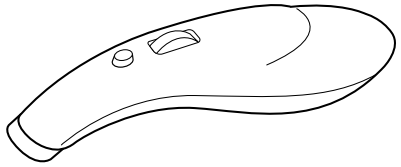
Seecode Thermo Insect Bite Device

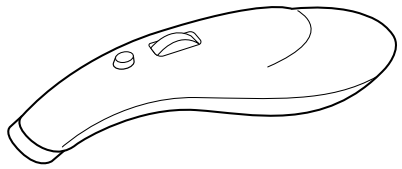EN I SEECODE Thermo Insect Bite Device Instruction manual

Thank you for purchasing our SEECODE Thermo Insect Bite Device, which helps you to quickly and easily prevent itching, pain and swelling caused by insect bites. To ensure that you can enjoy your product for a long time and make full use of it, please read these operating instructions carefully before using the SEECODE insect bite device to avoid damage due to incorrect operation.

Please pay special attention to the safety information. If you give the insect bite device to a third party, this instruction manual should also be handed over. Please keep the operating instructions for future

Never immerse the insect bite device in water.
Do not operate the device if it is damaged.
Only use the unit for the purpose described in these instructions.
This appliance can be used by children only if they are supervised or have been instructed in the safe use of the appliance and understand the hazards involved.

Do not use on or near mucous membranes, wounds, eyes or other sensitive body parts.
The device is not intended for use by persons (including children) who have little or no sensitivity to heat.
The insect bite device must not be used by persons with reduced

- and itching after an insect bite.
- The concentrated heat naturally destroys the proteins and thus eliminates the toxins that are the cause of the discomfort associated with the sting.
- For the safety of the users, this product is without chemical agents and additives. The material used is a material with thermal conductivity.

Safety instructions

This appliance is for household use only.

Technical data

Material ABS
Power supply 2 AA batteries (not included)
Generated heat approx. 40°C to 51°C

Instructions for use

First remove the protective film from the head section.

Using and/or replacing the batteries

- Open the battery compartment on the back of the unit.
- Insert 2 AA-Mignon batteries, observing the polarity (plus/minus).
- Close the battery compartment.

When the red LED lights up, this means that the batteries must be replaced.

sensation is different.
Always test before use on children.

- Increase the strength carefully and according to your sensitivity.
- After the itching has disappeared after approx. 2 to 5 seconds, put the device down and switch it off using the switch-off button.
- If necessary, the treatment can be repeated after 2 minutes.

The batteries must be inserted correctly. Please observe the polarity (plus/minus).
Incorrectly inserted batteries can damage the unit, they can leak and, in extreme cases, they can cause fire.
and in extreme cases it can lead to fire and explosion.
In case of malfunction, insert new batteries.

Always replace all batteries at the same time. Always use two of the same batteries (rechargeable, alkaline or other).
Remove used batteries directly from the unit.
We recommend that you remove the batteries if you are not going to use the appliance for a long time, as they may leak and damage your appliance.

Guarantee/Warranty

The warranty period is 24 months from the date of purchase. This includes a 6-month guarantee from the date of purchase. The territorial scope of the warranty protection is Europe-wide.

The proof of purchase is valid as proof of this. The product has been thoroughly checked for quality and functionality before dispatch. The guarantee is excluded in the following cases:
- in the event of modifications to the system without our approval
- in the case of use other than that stated above
- in the event of damage caused by falling

CE Declaration of Conformity

 Mobiset GmbH hereby declares that this product is in conformity with the essential requirements of the applicable European directives.

General information

Copyright

This document is protected by copyright. Any duplication or reprinting, including excerpts, as well as the reproduction of illustrations, even in modified form, is only permitted with the written consent of the manufacturer.

reference.

Our product is intended for adult use only. Please keep this product out of the reach of children for their safety.

WARNING

Children must not use the appliance without adult supervision.
Cleaning must also not be carried out by children.

Before using the insect bite device for the first time, remove the protective film from the head part.
The metal head heats up during use. Caution is advised.


physical, sensory or mental capabilities or by inexperienced or unformed persons. Unless they are supervised by persons responsible for their safety or have received instructions from them on how to use this appliance.

Important notes

- Only apply to skin areas that have insect bites. If the itching is still slightly noticeable afterwards, please repeat the treatment after 2 minutes at the earliest, but not more than 5 times an hour.
- This device uses heat to quickly and permanently relieve pain

Use the appliance only in accordance with the instructions in this manual.

Do not attempt to repair the appliance yourself.
It must be repaired by a customer service centre or a qualified repair centre.
Do not place or use the appliance near heat sources to avoid deformation of the plastic parts.

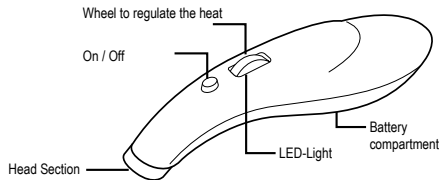Application

- Switch on the unit with the power button (clicks noticeably), the blue LED light flashes, indicating that the unit is working correctly.
- Place the head of the insect bite device on the insect bite.
- Use the dial to adjust the strength of the heat (from 1-low to 5-high).

Recommendation: Always start with strength 1.
Strength for children: 1 to 2 and for adults 2 to 5.
Attention. Depending on the person and body area, the temperature

Cleaning and care

Clean the appliance only with a dry or slightly damp cloth that does release any lint.
Do not use aggressive cleaning agents to clean the unit.

Handling batteries

Please take the used batteries to your local dealer. Do not dispose of them with normal household waste.
Limit the use of batteries, as disposal pollutes the environment.

Do not throw them into a fire, risk of explosion.
Charging the batteries and replacing batteries may only be carried out by adults.

When making a claim under the guarantee/warranty, please return the complete product in transport-safe packaging to the point of sale from which the product was purchased. Describe the detected defect in the product and enclose the original proof of purchase with the shipment. The seller reserves the right to remedy the defect before refunding the purchase price. Necessary transport costs will be reimbursed to the customer if a justified claim under the statutory guarantee/warranty is established. The use of the legal rights for defects is free of charge. These are not limited by the guarantee.

Notes on environmental protection

 The packaging materials used can be recycled. Dispose of packaging materials that are no longer needed in accordance with the locally applicable regulations. This product must not be disposed of in normal household waste within the European Union. Dispose of the appliance in accordance with regional regulations.

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