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# EN | Seecode Foot massager with EMS technology

## Instruction manual

Thank you for purchasing our SEECODE Foot Massager with EMS Technology, which allows you to send controlled electrical impulses directly to the muscles. Repeated contractions can relax muscles, reduce inflammation, prevent muscle atrophy, accelerate muscle healing and stimulate muscle growth.

In order to enjoy and make full use of your product for a long time, please read this instruction manual carefully before using the SEECODE foot massager with EMS technology to avoid damage due to incorrect operation. Please pay special attention to the safety information. If you give the SEECODE foot massager with EMS technology to a third party, this instruction manual should also be handed over. Please keep the operating instructions for future reference. Our product is intended for adult use only. Please keep this product out of the reach of children for their safety.

## WARNING

Children must not use the appliance without adult supervision.

Cleaning must also not be carried out by children.

Safety instructions


This appliance is intended for household use only. Only use the appliance in accordance with the instructions in this manual. Do not attempt to repair the appliance yourself. It must be repaired by a

It must be repaired by a customer service centre or a qualified repair centre.

## Scope of delivery

1. Massage mat
2. USB power cable
3. Small spray bottle to moisten the mat
4. Manual
5. Storage tube

## Description of the product

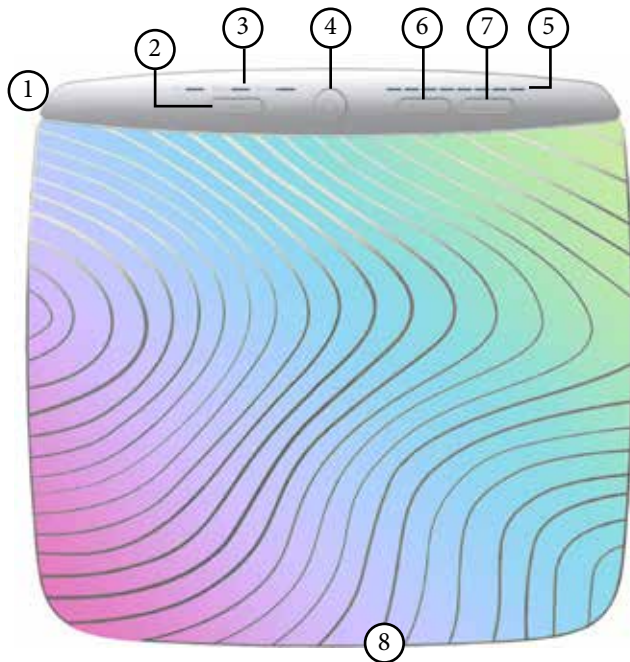1. Control unit
2. Button for setting the modes
3. Mode indicator
4. On / Off switch
5. Strength indicator
6. Strength minus
7. Strength plus
8. Massage mat

## Instructions for use

First remove the foil from the back of the mat.

The SEECODE foot massager with EMS technology is an electrical muscle stimulation device. It is equipped with 3 Hertz frequencies (HZ) from 30Hz to 10000Hz and

also 8 intensity levels that can be adjusted according to your comfort. The Hertz frequency is decisive for the EMS application. When selecting the frequency, be aware that there are individual differences.

0 Low Hertz frequencies mainly appeal to the slowly reacting red muscle fibres. For strength or endurance training, you should choose muscle stimulation in this frequency range.

0 Higher Hertz frequencies stimulate the fast contracting white muscle fibres. The muscle is deliberately overstrained and is thus stimulated to build up.

1. Open the massage mat (the standing surface for the feet), spread it out in front of you and stand on the mat so that the bare skin of the feet or other body parts touch the electrodes (curved lines in the mat) of the massage mat. Please do not wear socks or shoes when using on your feet.
2. Three hertz frequency modes are available, the mode is selectable from 30Hz to 10000Hz. The green indicator light shows the selected mode.

3. Now set the intensity of the impulses. You can choose between 8 levels (level 1 to level 8), depending on the level selected, the indicator lights in the control unit light up.

weak (◀)

strong (▶)

Always start with a low setting and increase it as you feel comfortable.

## Notes on use

1. The foot massager is suitable for use on the feet, calves, thighs, hips, arms and palms.
2. Place the product on a solid surface during use.
3. The usage time is 5 minutes, the application stops automatically after 5 minutes. The foot massager must be restarted for another application.
4. The silver lines on the massage mat are the current conducting part of the EMS, the skin must have direct contact with it.
5. If necessary, you can use the small spray bottle to wet the floor mat with a little water to increase the intensity of the EMS impulses.
6. Folding or creasing the mat may damage the product. Always roll up the mat carefully. Use the storage tube for transportation or storage.
7. Please fully charge the product before use.

## Technical data

Input voltage: 5V/0.5A

Rated power: 1W

Power supply : rechargeable lithium battery

Connection type : Type-C

Main material: ABS and TPU

Product weight: 220g

Product size: 383 x 357mm

Battery capacity: 300mAh

## Care and storage

1. Clean the massage mat after use with a damp cloth, wipe off any remaining water stains and store it rolled in the tube only.
2. Please store it out of the reach of children and pets.
3. If you do not use the foot massager for a long time, please charge it fully before storing it.
4. Store it in a dry place at about 10° to 40°, avoid direct sunlight, high temperature and humidity.

## Precautions

**The foot massager is not suitable for the following people:**

- Patients with cerebral haemorrhage
- Patients with thromboembolism
- Women in pregnancy and during menstruation
- Patients suffering from tuberculosis
- Patients with pacemakers
- Patients with severe renal failure
- Patients with severe heart disease
- Patients with malignant tumours
- Persons with acute varicose veins
- Persons with dermatitis and skin infections
- Persons with osteoporosis
- Persons with severe diabetes

## Safety instructions

1. Keep young children away from the foot massager.
2. Do not use the appliance together with other therapeutic appliances, beauty appliances or electrical appliances.
3. Please do not use it if you notice any malfunction or feel any abnormality in your body. Make sure that no metal objects touch the massage mat while it is switched on.
4. Do not use one application for more than 5 minutes at a time, the usage time in a day should not exceed 30 minutes. Avoid continuous use of the same part of the body, as long-term and continuous use may cause muscle damage.
5. Do not use the appliance in places with high humidity, such as bathrooms.
6. Do not use it on damp or heavily creamed skin.

## Battery and charging

The foot massager has a built-in rechargeable battery. Please do not place it in a fire or on a radiator. Charge it regularly, but please do not use it in very low temperatures. This may cause damage. When the power of the foot massager is too low, the indicator light turns red and the unit needs to be charged.

During charging, the eight intensity level indicator lights flash and light up cyclically. When the flashing stops, they light up continuously, indicating that the charging process is complete.

The foot massager cannot be used during the charging process.

Please only use a suitable charger for charging.

After about 2 hours of full charging, the total operating time is up to 8 hours.

Depending on the intensity of use and the environment in which you use the foot massager, the usage time may vary.

Please fully charge the foot massager when you use it or do not use it for a long time.

If the USB cable is damaged or the charging connector is loose, please stop using the foot massager.

### **Troubleshooting**

If a fault other than those listed below occurs and the foot massager cannot be used, never attempt to disassemble or repair the unit yourself to avoid danger.

### **Troubleshooting**

#### **The foot massager does not start**

Check if the foot massager has power, if there is no power, please charge it fully first.

#### **The battery does not charge**

Clean the charging port carefully before charging without getting moisture on the port.

#### **The foot massager stops during use**

Check the usage time, it will stop automatically after 5 minutes, just restart.

#### **No EMS current noticeable**

The device must be in direct contact with the skin during use, please also keep the skin and massage mat clean and free of grease at all times.

#### **The EMS current is only slightly perceptible.**

Wipe the used area with a damp cloth.

#### **The EMS current is too strong**

The foot massager is set too strong. Always start with the lowest levels and then set them higher according to your comfort.

## Cleaning and care

Clean the device only with a dry or slightly damp cloth that does not shed lint. Do not use aggressive cleaning agents to clean the appliance.

## Warranty / Guarantee

The warranty period is 24 months from the date of purchase. This includes a 6-month guarantee from the date of purchase. We grant a 6-month guarantee on the battery.

The geographical scope of the warranty protection is Europe-wide. The proof of purchase is valid as proof of this. The product has been thoroughly checked for quality and functionality before dispatch. The guarantee is excluded in the following cases:

- in the event of modifications to the system without our approval
- in the case of use other than that stated above
- in the event of damage caused by falling

When making a claim under the guarantee/warranty, please return the complete product in transport-safe packaging to the point of sale from which the product was purchased. Describe the detected defect in the product and enclose the original proof of purchase with the shipment. The seller reserves the right to remedy the defect before refunding the purchase price. Necessary transport costs will be reimbursed to the customer if a justified claim under the statutory guarantee/warranty is established. The use of the legal rights for defects is free of charge.

## CE Declaration of Conformity

Mobiset GmbH hereby declares that this product is in conformity with the essential requirements of the applicable European directives.

## General information

### Copyright

This document is protected by copyright. Any duplication or reprinting, including excerpts, as well as the reproduction of illustrations, even in modified form, is only permitted with the written consent of the manufacturer.

### Notes on environmental protection

The packaging materials used can be recycled. Dispose of packaging materials that are no longer needed in accordance with the locally applicable regulations.

This product must not be disposed of in normal household waste within the European Union. Dispose of the appliance in accordance with regional regulations.

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Mobiset GmbH - Rösrather Straße 333 - 51107 Köln