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# EN | SEECODE Mini Massage Gun

## Instruction manual

Thank you for purchasing our SEECODE Mini-Massage Gun, which can be used to quickly and pleasantly relieve pain and tension in muscles, sprains and strains as well as swelling. Thickened connective tissue and fascia can be loosened, muscle fatigue eliminated. To ensure that you enjoy your product for a long time and can use it to the full, please read these operating instructions carefully before using the SEECODE mini massage gun to avoid damage due to incorrect operation. Please pay special attention to the safety information. If you give the mini massage gun to a third party, this instruction manual should also be handed over. Please keep the operating instructions for future reference. Our product is intended for adult use only. Please keep this product out of the reach of children for their safety.

## WARNING

Children must not use the appliance without adult supervision.

Cleaning must also not be carried out by children.

## Safety instructions


This mini massage gun is intended for household use only. Only use the appliance in accordance with the instructions in this manual. Please do not attempt to repair the appliance yourself. It must be repaired by a customer service centre or a qualified repair centre.

## Scope of delivery

1. Massage gun with integrated battery
2. 4 massage heads
3. Charging cable Type-C / USB
4. User manual


## Product overview
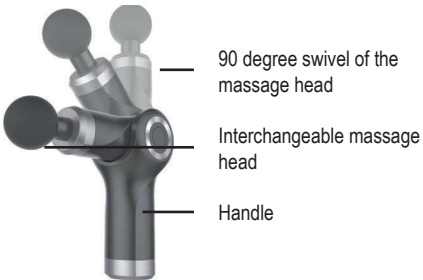
## Technical data

Battery:	2000mAH high capacity lithium battery
Running time:	up to 8 hours, depending on massage speed
Charging time:	approx. 3.5 hrs.
Interface type:	Type-C interface
Speed settings:	1600 / 2100 / 2600 / 3200 RPM
Noise level:	max. 50 dB
Number of massage heads:	4 pieces, cylindrical, round, flat, U-shaped
Housing material:	ABS plastic
Net weight:	0.41 kg
Product size:	10.8 x 4.3 x 15.5cm

## Massage heads
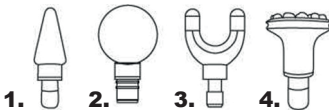
### 1. Conical head

For trigger points, pressure-sensitive hardenings in the tissue, precise point massage

## 2. Spherical head

For massaging large muscle groups such as arms and upper, lower legs


## 3. U-shaped head

For the neck area, spine and the Achilles tendons

## 4. Flat knobbed head

For all parts of the body, deep tissue massage

### Putting on and changing the massage heads

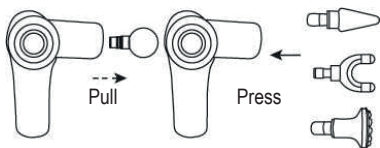
Select a massage head and place it on the appliance. To replace, pull the massage head you want to replace with a slight twisting motion and put on another one.

### Adjusting the massage head

### Indicator light (blue) and signal light (green and red)

Massage levels

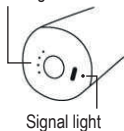The indicator light shows the massage levels. As the massage intensity increases, one more lamp lights up on the underside of the unit.


lights up on the underside of the unit.

When the unit is switched on, the signal lamp flashes the signal lamp flashes green once if the power is sufficient. The signal lamp lights up green 3 times every 5 seconds when the unit is fully charged. If the battery is less than 60% charged, it flashes 2 times every 5 seconds. If the battery is less than 30% charged, it flashes only once every 5 seconds.

The signal lamp lights up red when the power is less than than 10%.

The signal lamp lights up red when the battery is being charged.

## Charge

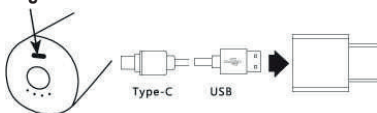*Insert the USB plug into the power adapter (not included)*

Please ensure that the device is fully charged before using it for the first time. Connect it to a suitable charger (DC5.0 ± 0.25V) using the USB cable provided.

The signal lamp lights up red when charging and green when fully charged (approx. 3.5 hrs.). Disconnect the charger from the massager when charging is complete to prolong the life of the battery.

If the red signal light comes on during use, recharge the massage gun as soon as possible.

The massage gun cannot be used while charging.


## Application

Switching on / off

Press and hold the start button for three seconds. The signal lamp and the indicator lights light up and the unit starts. To switch off, press and hold the start button again for three seconds until the unit switches off.

## Setting the massage speed

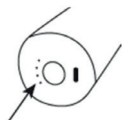After switching on, the first level is automatically switched on. Press the start button again to switch up one level at a time. The highest level is the fourth level.

Press the start button once in the fourth level to return

to the first level.

**Sequence:**

Switch on = 1st stage

Press start button = 2nd stage

Press start button = 3rd stage

Press start button = 4th stage

Press start button = 1st stage etc.

Safety note: To protect bones and muscles, do not massage the same area for too long. In any case, the appliance switches off automatically after 15 minutes.

**Important health note**

Do not use this product in the following cases/areas:

- Women during pregnancy
- Persons with heart disease, e.g. with a pace-maker
- People with cancer
- People with blood stasis
- People with tumours

If you have had an operation within 90 days, the product must not be used on this operative site of the body and in the vicinity of approx. 10 cm

Caution for sensitive areas: Head, face, cervical vertebrae and vertebral bodies of the spine.

Avoid positions close to the bone with little muscle coverage: e.g. the back of the foot and the back of the hand.

**Cleaning**

Always unplug the unit from the mains before cleaning. Never operate or charge the device unattended. Do not place or drop in water or other liquids.

Clean the unit only with a dry or slightly damp cloth that does not shed lint. Do not use aggressive cleaning agents to clean the unit!

**Disposal of the unit**

The European Directive on waste electrical and electronic equipment (WEEE) stipulates that this equipment must not be disposed of with household waste. Used equipment must be disposed of via separate collection to increase the reuse of the various components, thus protecting health and the environment.

Please dispose of the rechargeable batteries in the battery container or take them to a specialist shop. Never dispose of batteries in normal household waste.

### **Guarantee/Warranty**

The warranty period is 24 months from the date of purchase. This includes a 6-month guarantee from the date of purchase. We grant a 6-month guarantee on the rechargeable battery.

The territorial scope of the warranty protection is Europe-wide. The proof of purchase is valid as proof of this. The product has been thoroughly checked for quality and functionality before dispatch. The guarantee is excluded in the following cases:

- in the event of modifications to the system without our approval
- in the case of use other than that stated above
- in the event of damage caused by falling

When making a claim under the guarantee/warranty, please return the complete product in transport-safe packaging to the point of sale from which the product was purchased. Describe the detected defect in the product and enclose the original proof of purchase with the shipment. The seller reserves the right to remedy the defect before refunding the purchase price. Necessary transport costs will be reimbursed to the customer if a justified claim under the statutory guarantee/warranty is established. The use of the legal rights for defects is free of charge.

### **CE Declaration of Conformity**

Mobiset GmbH hereby declares that this product is in conformity with the essential requirements of the applicable European directives.

### **General information**

#### **Copyright**

This document is protected by copyright. Any duplication or reprinting, including excerpts, as well as the reproduction of illustrations, even in modified form, is only permitted with the written consent of the manufacturer.

## Notes on environmental protection

The packaging materials used can be recycled. Dispose of packaging materials that are no longer needed in accordance with the locally applicable regulations.

This product must not be disposed of in normal household waste within the European Union.

Dispose of the appliance in accordance with regional regulations.

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