Waden- und Oberschenkel-Massagegerät
Calf and thigh massager
Masseur de mollets et de cuisses
Kuit- en dij massagetoestel

DE (2) | EN (12) | FR (22) | NL (32)

EAN 4026924204221

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EN I SEECODE calf and thigh massager

Thank you for purchasing our SEECODE calf and thigh massager, which you can use to massage your calves and thighs separately. The ergonomic design adapts well to the calves and thighs and therefore also supports the compression massage.

The calf and thigh massager massages the muscles by alternating the pressure exerted by the dynamic air cushions. This compression massage promotes blood circulation in general and supports the leg muscles - even after a workout, for example.

To ensure that you can enjoy your SEECODE calf and thigh massager for a long time and use it to its full potential, please read these operating instructions carefully and pay particular attention to the safety information before using the SEECODE calf and thigh massager to avoid damage caused by incorrect operation.

If you give the SEECODE calf and thigh massager to a third party, these operating instructions should also be handed over. Please also keep the operating instructions for future reference.

Important safety instructions for your health, please observe!

Disclaimer - Health risks when using the massager

- Use of the calf and thigh massager is at your own risk. The device is not intended as a substitute for medical treatment or physiotherapy. In the event of existing health complaints, particularly in the area of the muscles, nerves or blood circulation, a doctor or medical professional should be consulted before use.
- If there is a suspicion of thrombosis, high blood pressure, diabetes or heart disease, the doctor treating you must also be consulted before use.
- The manufacturer/distributor accepts no liability for damage to health caused by improper use or failure to observe the instructions for use. Persons with acute or chronic illnesses, pregnant women and people with medical implants (e.g. pacemakers) should only use the device after consulting a doctor.
- This product is not suitable for patients who have recently undergone surgery or who suffer from burns, skin inflammation, sensitive skin or purulent diseases.

Technical data

Maximum leg circumference: with wide cuff circumference for the thighs up to 74 cm,
with narrow cuff circumference for the calves up to 65 cm

Batteries: 1 x 2500 mAh lithium-ion battery each (integrated battery)

Remote control battery: 1 x CR2032 each (replaceable)

Charging cable: USB-C / USB-A connection cable (Power supply unit 5V, 2A not included)

Power: 8 watts

Input voltage: 5 volts

Massage intensity: 3 massage modes, 3 heat modes and 5 intensity modes adjustable

Heating levels: 37°C, 42°C and 47°C

Automatic switch-off: after 20 minutes

Size: 70 x 30 cm with narrow connecting piece to the massage unit for the calves and
79 x 30 cm with wide connecting piece to the massage unit for the thighs

Weight: 900 grams (both massage units)

Outer material: polyester

Charging time: approx. 3 hours when the battery is fully discharged

Running time: approx. 1.5 hours, depending on the selected mode

Scope of delivery

2 x calf and thigh massagers with permanently mounted base unit

One remote control per device, can also be used on the base unit as a control unit

2 connectors for use on the calves or thighs

Battery charging cable (Power supply unit not included)

Operating Instruction


Important instructions, please read and observe!

- The massagers are labelled with an L for use on the left leg and an R for use on the right leg. Please use the version with the L for the left leg and the one labelled R for the right leg.
- Keep the calf and thigh massager out of the reach of children. Do not allow children to play with this product as a toy.
- Do not rinse the device under running water and do not place it in water.
- Do not use the neck and shoulder massager if it is defective in any way.
- Do not use the product during sleep, after extensive fluid intake or after drinking alcohol. If excessive pain or skin discomfort occurs during use, stop using the massager immediately and consult a doctor.
- Stop using the product immediately if it becomes deformed or overheated.
- Do not use the product in a damp environment, near fire or if there are wounds on the skin.
- Even at the lowest temperature setting, people with very sensitive skin may suffer burns. When using the heating function, make sure that there are no skin abnormalities, especially before extending the usage time. As individual skin conditions vary, the perceived

temperature may differ.

- Disconnect the mains plug after the device has been charged.
- Please do not stand on the device.
- If the device is damaged, stop using it immediately to avoid accidental damage.
- Do not use a damaged power cable or plug. If the cable or plug of the power supply unit is damaged, have it replaced by a specialist.

Charging the battery

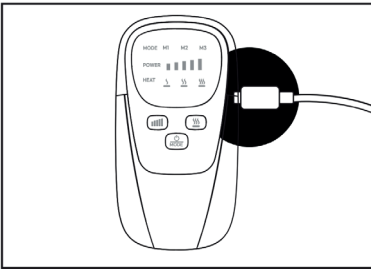1. The charging connection is located on the side of the respective permanently installed base unit.

Only use the enclosed charging cable for charging.
(Power supply unit not included)

2. The product cannot be used during the charging process.

During the charging process, you can check the charging status on the base unit. While the battery is charging, the symbols for 25%, 50%, 70% and 100% charge light up.

3. depending on the charge status of the battery, it may take up to 3 hours until it is fully charged.

Important note: Basic switching on and off of the massager

Remove the remote control from the permanently mounted base unit.

Always switch on the massager before using the remote control by pressing and holding the ON/OFF button on the base unit for a few seconds.

When the device is switched on, the standard program with the settings starts immediately:

Mode 1

Massage intensity 1


Heat level 1

To switch off the massager, press the ON/OFF button on the base unit again for a few seconds.

After switching off the device, the air slowly escapes from the air chambers. DO NOT try to press the cuff to expel the air!

Operating the calf and thigh massager

1. Putting on the massager

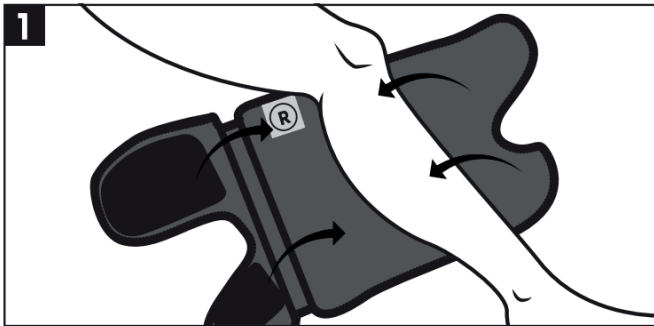First check which side you can put your massager on by following the instructions (L or R) on the product label.

Note: For use on the calf, please select the narrower connecting piece; for use on the thigh, please select the wider connecting piece. These are attached to the massage unit using a zip fastener.


Place the massager around your calf or thigh and adjust the fit to your leg circumference using the Velcro fasteners.

Note: Make sure the fit is comfortable and not too tight and position the massager so that the control unit is facing outwards.

Adjust the size with the Velcro fasteners according to the desired intensity level.

If the pressure from the massager becomes too strong for you, stop the application immediately by pressing the on/off button or immediately open the zip for the connecting pieces.

Operating the base unit

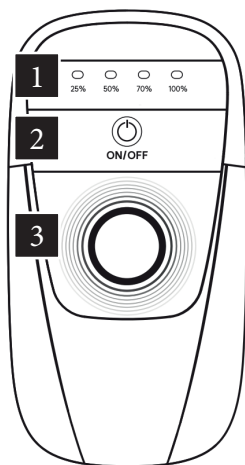
1. Charge level indicator from 25% to 100% charge

2. On/off button

This confirms with a “beep” that the massager is switched on. Standard settings are: Massage mode M1, intensity L1, heating off.

3. Remote control aiming point: When using the remote control, please point it at the base unit.


2. Switching on the remote control

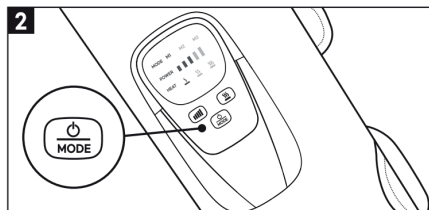Press the on/off/MODE button on the remote control for a few seconds.

The currently set mode settings light up in the display.

3. Making individual settings

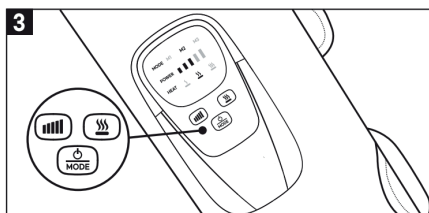
To set the mode, intensity and heating, use the corresponding buttons on the remote control.

- Setting the intensity: Use the intensity button to select the intensity in levels 1 to 5.
- Setting the heating: Use the heating button to select levels 1-3 or OFF. Note: The heating function cannot be used without using other functions.


- Selecting the preset massage modes: After switching on the device, use the mode button to switch between mode 1, mode 2 and mode 3.
Select the massage mode that is most comfortable for you.

4. Automatic switch-off

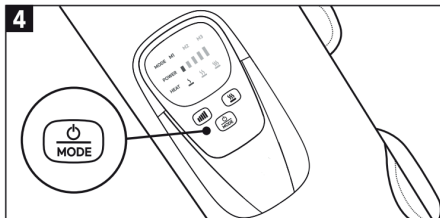The device switches off automatically after 20 minutes of operation.
To force a switch-off, press and hold the On/Off/ MODE button on the remote control for a few seconds.

5. Operating the permanently installed base unit

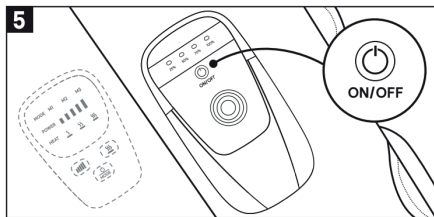You can only switch the massager on and off completely using the ON/OFF button on the base unit.

The default settings for the switch-on button on the base unit are: M1 mode, intensity level 1 and heating off.

To change the mode, set the intensity or control the heating, you must use the remote control.

Operation with the remote control

The remote control can be easily lifted off the base unit.

However, it can also remain on the base unit when using it for individual settings.

1. Massage mode display

Mode 1, Mode 2, Mode 3

2. Intensity display

Intensity 1, Intensity 2, Intensity 3, Intensity 4, Intensity 5

3. Heating intensity display


Intensity 1, Intensity 2, Intensity 3

4. Button for setting the intensity level




5. On/Off/Mode button for the default settings:

Mode 1, intensity level 3, heating level 1

6. Button for setting the heating intensity

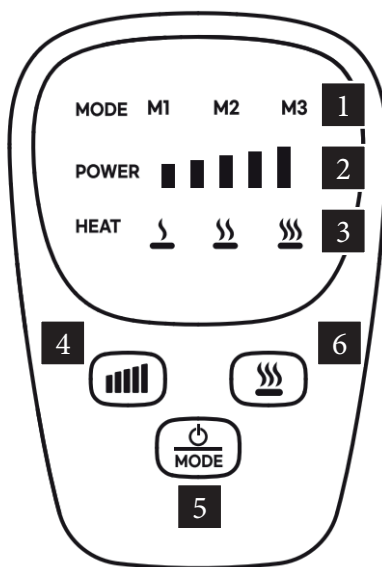Removing the remote control cover and inserting a battery (CR2032)

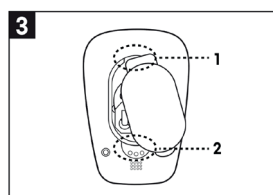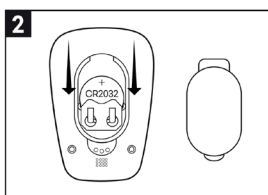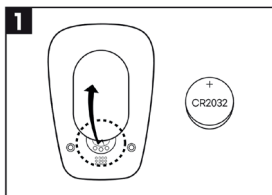1. Use the cut-out to open the back of the remote control and remove the cover. We recommend using a tool to remove the cover, taking care with your hands. Remove the battery.

2. Insert a new battery (CR2032) as shown in the illustration and then push it in from top to bottom. Please note the correct polarity (+ / -) of the battery.

3. Insert the wider tab on the cover into the opening (1) and then press on the lower part of the cover until you hear a click in the lower opening (2).

Cleaning and care

Remote control / main control unit

Wipe the remote control and the main control unit carefully with a soft cloth and store them in a safe and clean place.

Outer fabric

Never wash the calf and thigh massager in water.

Only wipe the soiled areas lightly with a damp cloth and a mild detergent if necessary. Do not use abrasive cleaning agents or harsh chemicals. Do not use any metallic cleaning objects.

General safety instructions for non-removable batteries

- Batteries may only be charged by children under adult supervision. Please keep away from children at all times.
- The electrolytes and electrolyte vapors contained in a rechargeable battery are harmful to health. Avoid direct contact at all times. If electrolytes come into contact with skin, eyes or other parts of the body, rinse immediately with plenty of fresh water. Then consult a doctor.
- Stop use or the charging process immediately if the battery swells, leaks, has a high temperature or emits smoke and strong odors.
- If you do not use the device for a longer period of time, charge it in between. Connected batteries can become deeply discharged and must then no longer be used or charged.
- Do not immerse batteries in water.
- Do not short-circuit the battery (connect + and - polarity).
- Do not disassemble or open batteries, as this can lead to burns from the electrolytes.
- Do not expose batteries to excessive heat, e.g. direct sunlight, fire or similar. Do not place them on a heater.
- Do not dispose of rechargeable batteries by destroying or burning them and do not dispose of them in household waste. Dispose of them in accordance with local regulations.
- Please note that the end user is legally obliged to return used batteries and that used batteries can be returned to dealers free of charge.

Guarantee/warranty

The warranty period is 24 months from the date of purchase. This includes a 6-month guarantee from the date of purchase. We grant a 6-month guarantee on the battery.

The geographical scope of the warranty cover is Europe-wide. The proof of purchase serves as proof. The product has been thoroughly checked for quality and functionality before dispatch.

The guarantee is excluded in the following cases

- if the product has been modified
- if it has been used in a manner other than that specified above
- if it has been dropped

If you make a claim under the guarantee/warranty, please return the complete product in transport-safe packaging to the point of sale from which the product was purchased. Describe the recognized defect in the product and enclose the original proof of purchase with the consignment. The seller reserves the right to subsequent fulfilment before reimbursement of the purchase price. Necessary transport costs will be reimbursed to the customer if a justified claim under the statutory guarantee/warranty is established. The utilization of statutory warranty rights is free of charge.

CE Declaration of Conformity

Mobiset GmbH hereby declares that this product is in conformity with the essential requirements of the applicable European directives.

General information

Copyright

This document is protected by copyright. Any duplication or reprinting, even in part, as well as the reproduction of the illustrations, even in modified form, is only permitted with the written consent of the manufacturer.

Notes on environmental protection

The packaging materials used are recyclable. Dispose of packaging materials that are no longer required in accordance with local regulations.

This product must not be disposed of with normal household waste within the European Union. Dispose of the device in accordance with regional regulations.

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