Lendenmassage-Gerät
Lumbar Massage Device
Appareil de massage lombaire
Lendenmassageapparaat

DE (2) | EN (12) | FR (22) | NL (32)

EAN 4026924204436

DE | SEECODE Lendenmassage-Gerät

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EN I SEECODE Lumbar Massage Device

Thank you for choosing the Seecode Lumbar Massage Device.

This device offers you various functions for effective massage of the lumbar region at any time. It helps to relieve pressure and fatigue in this area.

The massage device is ergonomically adapted to the human body shape and is therefore suitable for all users. An integrated air cushion lifts the device individually to gently stretch the lumbar spine.

The stretching is combined with heat, low-frequency pulses (EMS) and vibration massage.

For discomfort in the lumbar region, we recommend daily use for 10 to 30 minutes to relax the muscles. If you have severe discomfort or existing medical conditions, please consult your doctor before use.

Please read this instruction manual completely and carefully before using the device. This will help you avoid damage caused by improper use. Pay particular attention to the safety information.

If you pass the device on to a third party, please include these instructions. Keep them in a safe place for future reference.

Safety instructions

The device is intended for use by adults only.

To avoid skin irritation or damage, do not use for excessive periods of time.

Scope of delivery:

Lumbar massage device with wired remote control

Plate for stabilising the surface

USB-A/USB-C charging cable (5V / 2A power supply required, not included)

Operating instructions

Technical data

Dimensions: 42.5 x 19.5 x 7.0 cm

Weight: 1,400 grams

Battery supply: DC7.4 V

Battery capacity: 2000 mAh

Charging time: approx. 2 hours when battery is empty

Running time: approx. 2 hours, depending on function and mode

Charging voltage: 5 V / 2 A

Nominal power: 15 W

Static operating current <30 uA

Switches off automatically after 10 minutes.

Three heat modes: low – medium – high.

Functions:

Height adjustment function (air cushion) with 15 modes

Low-frequency pulses (EMS) with 20 modes

Heat function

High-frequency vibration

Operation of the lumbar massage device

Remote control

- 1 Heating intensity display
- 2 Height adjustment display
- 3 Height adjustment lock display
- 4 Increase height adjustment
- 5 On/off button and height lock button
- 6 Reduce height adjustment
- 7 Heating intensity selection button
- 8 Vibration intensity display
- 9 EMS intensity display
- 10 EMS mode display
- 11 Increase EMS intensity
- 12 EMS mode selection
- 13 Reduce EMS intensity
- 14 Vibration function button
- 15 Charging indicator light

On/off button/lock button

Switching on: Press and hold button 5 for a few seconds until you hear a beep.

Switching off: Press and hold button 5 again for a few seconds. A beep will confirm the process.

Lock function: Press button 5 briefly to lock the height adjustment.

After switching on, you can activate or adjust the respective function by briefly pressing the other function buttons.

Note:

If the maximum level is reached for a function, the mode automatically returns to the lowest level the next time you switch, except for height adjustment.

Description of functions

Height adjustment +

- Function: Increases the height of the massage unit.
- Levels: Maximum 15 levels.
- Operation: Pressing the button briefly increases the level by 1.

Height adjustment –

- Function: Lowers the height of the massage unit.
- Levels: Maximum 15 levels.
- Operation: Pressing the button briefly decreases the level by 1.

EMS mode selection

- Function: Selects the pulse mode.
- Modes: A, B, C, D.
- Operation: Pressing the button briefly switches to the next mode.
- Automatic: After mode D, the setting automatically returns to mode A.

EMS intensity +

- Function: Increases the strength of the low-frequency pulses.
- Levels: 20 intensity levels in all modes.
- Operation: Pressing the button briefly increases the intensity by 1 level.

EMS intensity –

- Function: Decreases the strength of the low-frequency pulses.
- Levels: 20 intensity levels.
- Operation: Pressing the button briefly decreases the intensity by 1 level.

Important note: The EMS function only works when all 4 electrodes are in contact with the body at the same time. Otherwise, the function is immediately deactivated and must be reactivated using the remote control.

Please always start the application at the lowest level.

Heating intensity selection

- Function: Switches the heating function on/off and selects the heating level.
- Levels: 3 heating levels – the higher the level, the higher the temperature.
- Operation: Press once to switch on the heating; press again to change the level.

Vibration function

- Function: Switches the vibration on/off and selects the vibration mode.
- Modes: 3 vibration modes – increasing intensity with each level.
- Operation: Press once to switch on; press again to increase the level.

Using the stabilization plate

The stabilization plate included in the scope of delivery can be placed under the device to ensure additional stability.

Recommended use:

In soft environments, such as on soft beds or sofas, the stabilization plate improves the massage experience.

Not required:

Use is not necessary on hard surfaces (e.g. floors, firm mattresses).

Charging the battery

To charge, plug the USB-A connector of the charging cable into a power adapter (not included) with a power output of 5 V / 2 A.

Connect the USB-C plug to the charging port on the side of the lumbar massage device.

Fully charge the battery before using it for the first time.

During charging, the charging indicator on the remote control display flashes; when fully charged, the indicator lights up continuously.

When the battery is empty, a full charge takes approx. 2 hours.


If the device is not used for a longer period of time (weeks/months), charge it beforehand to prevent deep discharge and associated battery damage.

Using the lumbar massage device

Preparation

- Make sure that the device is fully charged.
- Lay it flat on a suitable surface.
- Use the stabilization plate on soft surfaces (e.g. bed, sofa).
- Use a DC 5 V / 2 A charger (not included) to charge the device.

Positioning

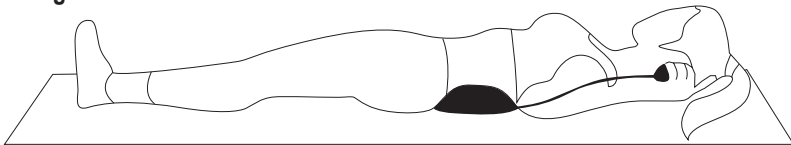- Lie down gently on the device.
- Adjust to your desired lying position.


Using the air cushion function

- Lift the massage device with the air cushion function to gently stretch the lumbar spine.

Activating functions

- Control the desired functions using the remote control.
- All functions can be combined with each other.
Note: The EMS function only works when all 4 electrodes are in contact with the body at the same time. Otherwise, it is immediately deactivated.
- Please always start the applications at the lowest level.

Ending the application

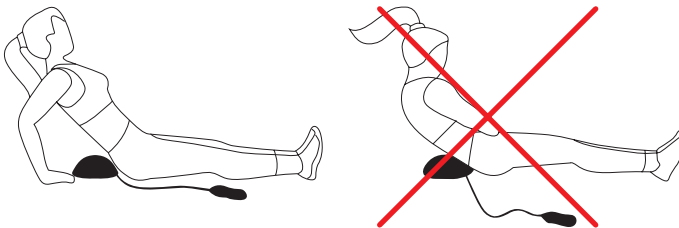- Deflate the airbag by setting the height adjustment to '0'.
- Remain lying down for 1–3 minutes after stretching.
- Then support yourself with your arms and stand up slowly and carefully.
- Caution: Avoid standing up abruptly after stretching to prevent straining your back muscles.

Important precautions, please note

The lumbar massage device must not be used:

- After lumbar spine surgery.
- By persons who are taking hormones over a long period of time and are extremely tired.
- Immediately after eating, drinking alcohol, bathing or intense physical activity.
- By pregnant women.
- By persons with implanted medical devices (e.g. pacemakers).
- By people with sensory disturbances (e.g. due to intense heat and/or electrical impulses on the surface).
- By people who are not sensitive to heat at all – in this case, special caution is required.

A doctor should be consulted before use in the following cases:

- Tumor diseases.
- Severe heart disease.
- High blood pressure.
- Severe skin diseases in the lower back area.
- Pronounced bone hyperplasia.
- Abnormal vertebral shapes or spinal curvature caused by accident/illness.
- Medically prescribed convalescence or unexplained physical complaints.
- High fever.
- Severe circulatory disorders in the lumbar region.

General warnings and safety instructions

Please stop the massage immediately if you experience unusual or severe pain during use – consult a doctor.

Stop use immediately in the event of device failure or unusual malfunctions.

Do not use in damp or dusty environments – risk of electric shock or device damage.

Protect from direct sunlight and do not place near heating appliances or ovens.

Do not use at ambient temperatures above 40 °C.

Do not use the device while it is charging.

Do not press any objects or fingers onto the airbag – risk of injury or product damage.

The device should always be firmly pressed against the skin during EMS application.

Do not use while sleeping.

The maximum recommended application time is 30 minutes per session.

Do not open the housing – risk of electric shock and damage to the device.

If the charging cable is damaged, be sure to use a suitable replacement.

The device is not suitable for use by persons (including children) with physical, sensory or mental impairments or lack of experience, unless under the supervision of a responsible person.

After use, unplug the power cord and allow the device to cool down before storing it in a plastic bag.

Store in a dry, cool place protected from moisture, high temperatures, direct sunlight and splashing water.

Troubleshooting

Possible Problems	Solution
The device does not charge when connected to the power supply	1. Check that the adapter is a plug-in power supply with 5V and min. 2A. 2. Check that the 5V socket is loose.
The device now heats up very slowly	The heating time of the device varies in different environments.
No EMS function noticeable	Ensure that the skin in the lumbar region is in direct contact with the 4 electrodes.
Discomfort when adjusting the height	Be sure to adjust the height set by air pressure to suit your individual needs. Some people's bodies react to even slight changes in height. This can be alleviated or eliminated by placing a pillow under the neck.

Cleaning

Only clean the device with a soft, slightly damp cloth. Do not wet the device or wipe it with a wet cloth – risk of water ingress.

Do not use abrasive cleaning agents or harsh chemicals. Do not clean with thinner, petrol or other volatile liquids.

Use a neutral cleaning agent for heavier soiling.

Safety instructions for handling permanently installed batteries

- Batteries may only be charged by children under adult supervision. Please always keep them away from children.
- The electrolytes and electrolyte vapors contained in a battery are harmful to health. Avoid direct contact at all costs. If electrolytes come into contact with skin, eyes or other parts of the body, rinse immediately with plenty of fresh water. Then consult a doctor.
- Stop using or charging the battery immediately if it swells, leaks, becomes hot, or emits smoke or a strong odor.
- If you do not use the device for a long period of time, charge it occasionally. Connected batteries can become deeply discharged and must then no longer be used or charged.
- Do not immerse batteries in water.

- Do not short-circuit the battery (connect + and - poles).
- Do not disassemble or open batteries, as this may cause burns from the electrolytes.
- Do not expose batteries to high heat, e.g. direct sunlight, fire or similar. Do not place them on a heater.
- Do not dispose of batteries by destroying or burning them, and do not dispose of them in household waste. Follow the locally applicable regulations when disposing of batteries.
- Please note the legal obligation of end users to return used batteries and the option of returning used batteries to retailers free of charge.

Warranty/Guarantee

The warranty period is 24 months from the date of purchase. This includes a 6-month warranty from the date of purchase. We provide a 6-month warranty on the battery.

The geographical scope of the warranty protection is Europe-wide. The proof of purchase serves as evidence of this. The product was thoroughly checked for quality and functionality before shipping.

The warranty is excluded in the following cases:

- if changes have been made to the system
- if the product has been used in a manner other than that specified above
- in the event of damage caused by dropping

To make a claim under the warranty/guarantee, please return the complete product in secure packaging to the point of sale where the product was purchased. Describe the defect you have identified in the product and enclose the original proof of purchase with the shipment.

The seller reserves the right to remedy the defect before refunding the purchase price. The customer will be reimbursed for any necessary transport costs if a justified claim under the statutory warranty/guarantee is established. The exercise of statutory rights in respect of defects is free of charge.

CE Declaration of Conformity

Mobiset GmbH hereby declares that this product complies with the essential requirements of the applicable European directives.

General information

Copyright

This document is protected by copyright. Any duplication or reprinting, even in part, as well as the reproduction of illustrations, even in modified form, is only permitted with the written consent of the manufacturer.

Environmental protection information

The packaging materials used are recyclable. Dispose of packaging materials that are no longer needed in accordance with local regulations.

This product must not be disposed of in normal household waste within the European Union. Dispose of the device in accordance with regional regulations.

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